



Ultrawave

hygea[®] 2850VM

HTM compliant ultrasonic cleaner



Features:

- Menu-driven LCD control
- Dual cycle traceability
(Electronic and/or hard copy)
- Electro-magnetic locking lid

Eastgate Business Park,
Wentloog Avenue,
Cardiff, CF3 2EY, UK

Tel: +44 (0) 845 330 4236
Fax: +44 (0) 845 330 4231
admin@ultrawave.co.uk
www.ultrawave.co.uk





Ultrawave

hygea[®] 2850VM

HTM compliant ultrasonic cleaner

Ultrawave's **Hygea 2850VM** ultrasonic cleaner meets the design considerations of HTM01-05 and HTM2030 for decontamination of dental instruments.

Designed specifically for use in the dental instrument decontamination process, the Hygea 2850VM is also commonly used in other healthcare sectors.


Combining high performance transducers and a 32 litre robust stainless steel tank, the **Hygea 2850VM** provides excellent, reliable performance for thorough cleaning of dental instruments.

The electro magnetic locking lid prevents instrument removal before being fully cleaned and increases safety of operators.



Every cleaning cycle can be recorded either in hard copy via the integral printer or electronically through a combination of the built in MC card port.

Supplied complete with a stainless steel basket and 2 x 3 tier basket tray system, the **Hygea 2850VM** can clean up to 90 standard dental instruments in a 10 minute cycle

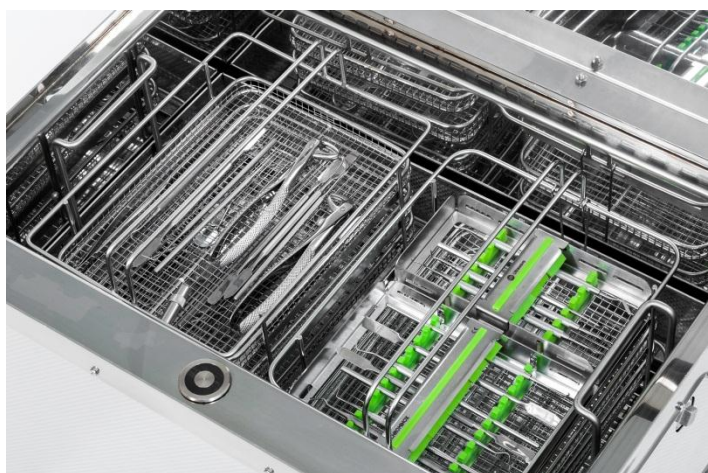

Eastgate Business Park,
Wentloog Avenue,
Cardiff, CF3 2EY, UK

Tel: +44 (0) 845 330 4236
Fax: +44 (0) 845 330 4231
admin@ultrawave.co.uk
www.ultrawave.co.uk





Ultrawave

hygea[®] 2850VM

HTM compliant ultrasonic cleaner

1. Rapid & effective instrument cleaning

Cleans up to 12 cannulated instrument per cycle

2. Menu-driven LCD control

Fully programmable split-level security gives excellent process control.

3. Digital cycle traceability

The MC card port is a brand new feature on the Hygea 2850VM and allows digital cycle traceability and validation of every cleaning cycle

4. Integral validation printer

Offers a hard copy print out of every cycle if required.

5. Electro-magnetic locking lid

Prevents instrument removal before being fully cleaned and increases safety of operators.

6. Basket shelf


Minimises wear & tear of the radiating surface

7. Stainless steel basket

Supplied with every machine.

8. 3 tier basket tray system

2 x systems are supplied with every machine.

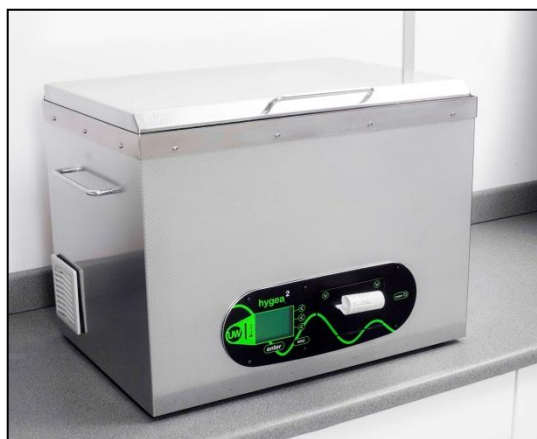

Eastgate Business Park,
Wentloog Avenue,
Cardiff, CF3 2EY, UK

Tel: +44 (0) 845 330 4236
Fax: +44 (0) 845 330 4231
admin@ultrawave.co.uk
www.ultrawave.co.uk





Ultrawave

hygea[®] 2850VM

HTM compliant ultrasonic cleaner

CONSTRUCTION CASE	Rigid stainless steel
CONSTRUCTION TANK	Stainless steel
OVERALL DIMENSIONS	520 x 330 x 330mm (WxDxH)
TANK DIMENSIONS	500 x 300 x 230mm (WxDxH)
BASKET DIMENSIONS	485 x 285 x 75mm (WxDxH)
CAPACITY	30 Litres
ULTRASONIC POWER	500W
OPERATING FREQUENCY	30-40KHz
HEATER POWER	750W
Drain Valve	3/4" BSP Valve
Control	Digital, menu driven LCD
Validation	<ul style="list-style-type: none">• Date & Time• Cycle no• Practice & Operator names• Set & achieved cycle parameters
Operating features	<ul style="list-style-type: none">• Electro-magnetic locking Lid prohibits process interruption• MC card data-logger for digital traceability• Integral printer for hard copy validation
MAINS CONNECTION	220v-240v (50-60Hz)

Eastgate Business Park,
Wentloog Avenue,
Cardiff, CF3 2EY, UK

Tel: +44 (0) 845 330 4236
Fax: +44 (0) 845 330 4231
admin@ultrawave.co.uk
www.ultrawave.co.uk

